


1 Preparation of the restoration*

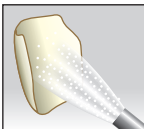


Sandblasting (Alumina oxide 50-100 µm)	
Adhesive surface	Pressure
Ceramics (silicate glass, lithium silicate glass), Composite resin	Approx. 0,1-0,2 MPa (Approx. 1-2 kgf/cm ²)
Ceramics (alumina, zirconia)	Approx. 0,2-0,3 MPa (Approx. 2-3 kgf/cm ²)


Roughen the restoration surface with an abrasive, or sandblast as described in Table, then rinse with water and dry.
Hydrofluoric acid etching is recommended for glass ceramic restoration.

*Translucent restoration with a thickness of less than 2mm.


2 Try-in

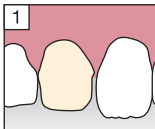

Remove the temporary sealing or provisional restoration from tooth, then rinse with water and dry.

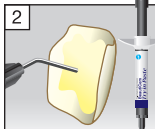


Select the proper shade of BeutiCem Try-in Paste and apply to the adhesive surface of the restoration.

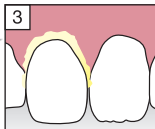


Place the restoration onto the tooth to confirm the shade.

Try-in

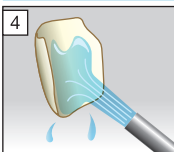


After try-in, thoroughly rinse off the Try-in Paste remaining on restoration surface.

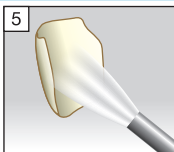


Dry the restoration.
In case of shade mismatch, retry with another shade of Try-in Paste from step 2-2.

3 Pretreatment of the restoration

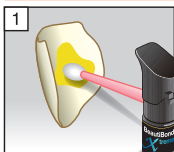

Apply BeutiBond Xtreme and air dry.
(Blow gentle air:3sec.→Blow more strongly).

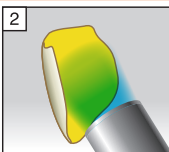

Light-cure.
LED : 5 sec,
Halogen : 10 sec.

※In the case where the adhesive surface is large, light-cure in segments.


4 Pretreatment of the tooth




Rinse with water.




Dry the tooth and confirm that there is no residual Try-in Paste on the tooth surface.

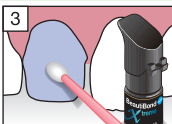

Apply BeutiBond Xtreme.

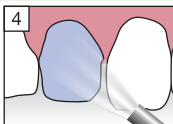

Air dry.
(Blow gentle air:3sec.→Blow more strongly).

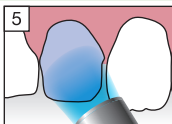


Light-cure,
LED : 5 sec,
Halogen : 10 sec.


※In the case where the adhesive surface is large, light-cure in segments.

Etch the tooth surface if necessary and then rinse thoroughly with water and air dry.


5 Placement/Seating of the indirect restoration

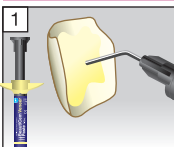


Apply BeutiCem Veneer to the adhesive surface of restoration.

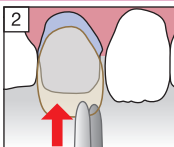

Place it on the tooth and apply gentle pressure to extrude excess cement from the margins.

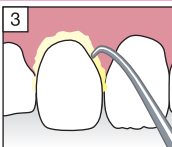


Tack-cure the paste for 2-3 sec and remove the excess cement.

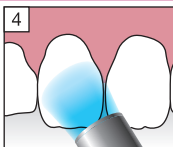

Light-cure.
LED : 10 sec,
Halogen : 20 sec.

※In the case where the adhesive surface is large, light-cure in segments.

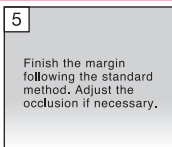

Finish the margin following the standard method. Adjust the occlusion if necessary.