


# Systematic Contouring, Finishing and Polishing of Composites

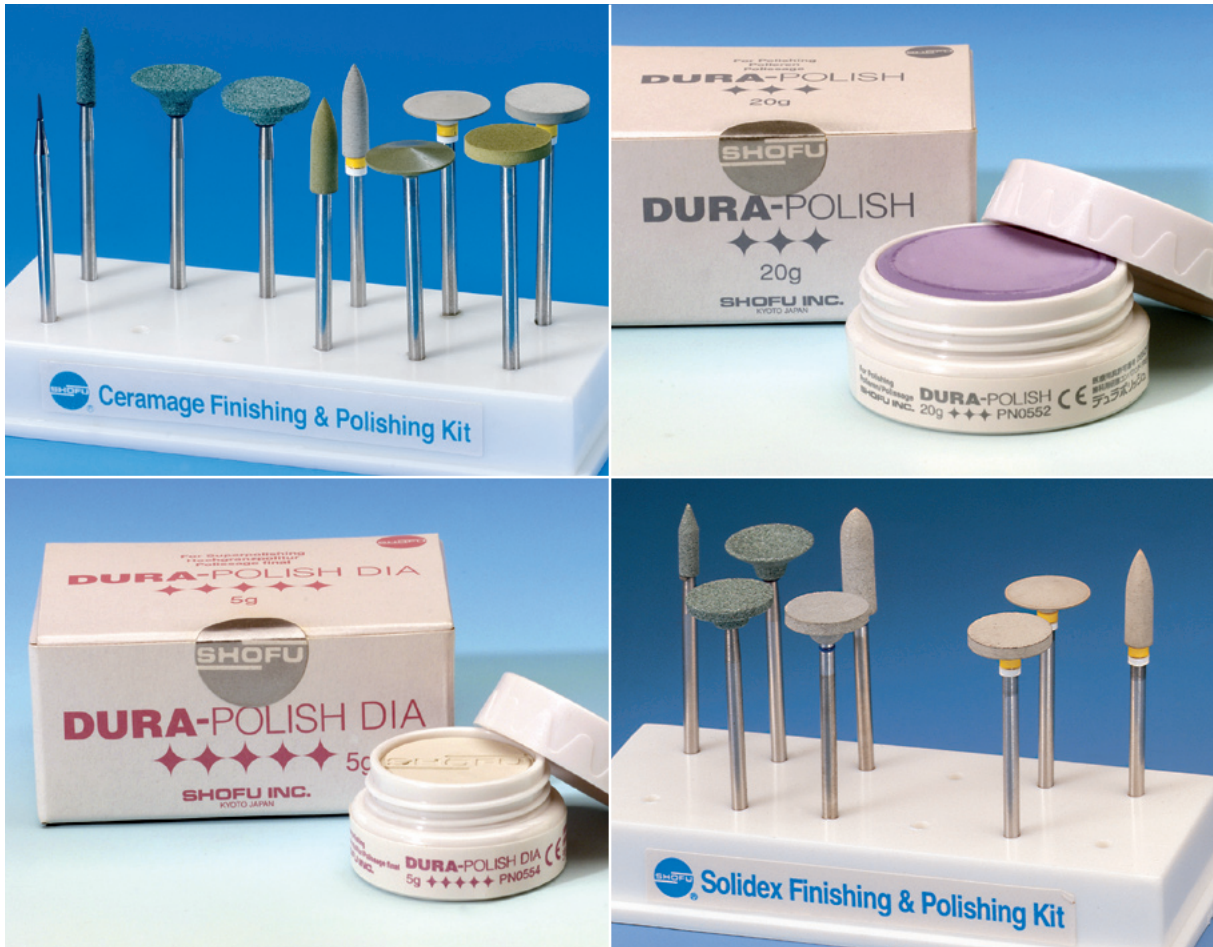

## Time-Saving Instrument Kits and Polishing Pastes

SHOFU's light-cured restoratives and ceramic-reinforced crown & bridge composites such as the micro-hybrid composite systems CERAMAGE and SOLIDEX offer excellent polishability for high and long-lasting surface quality and resistance to plaque build-up.

This is why systematic contouring, finishing and polishing of crowns / bridges, inlays / onlays or veneers after polymerisation are so important.

Our well-matched rotary instruments and highly filled polishing pastes allow you to easily, quickly and consistently create high-quality surfaces, even in hard-to-reach areas.



Proven Products for Better Dentistry

# SOLIDEX and CERAMAGE

## Finishing & Polishing Kits

The well-matched components of these two instrument kits allow you to systematically finish and pre-polish ceramic-reinforced micro-hybrid composites.



### Solidex Finishing & Polishing Kit

- 3 Dura-Green Stones
- 2 SoftCut PB Pre-polishers
- 3 CompoMaster Polishers

PN 0328



### Ceramage Finishing & Polishing Kit

- 1 Robot Carbide Fissure Bur
- 3 Dura-Green Stones
- 3 CompoMaster Coarse Pre-polishers
- 3 CompoMaster Polishers

PN 0333



### Robot Carbide Fissure Bur

The tip of this bur has the shape of a triangular pyramid which makes it particularly suitable for pits and fissures. Composite surfaces are gently adjusted and smoothed in one step with only minimal heat generation.



**Dura-Green** silicon carbide stones are used to adjust and contour tooth shape and surface.



**SoftCut PB** silicone pre-polishers are ideal for pre-polishing rough composite surfaces and pre-smoothing fine anatomical details.



**CompoMaster Coarse** diamond-impregnated silicone polishers contour and smooth rough composite surfaces in one step, depending on the contact pressure applied.



**CompoMaster** silicone polishers for finishing and polishing are impregnated with ultrafine diamond particles (approx. 6 µm) specially matched to the particle size of the ceramic fillers used in CERAMAGE and SOLIDEX.

## Dura-Polish & Dura-Polish DIA High-Performance Polishing Pastes



**Dura-Polish**, impregnated with aluminium oxide particles (73 % by weight), is used for pre-polishing and creates extremely homogeneous surfaces, even in hard-to-reach areas. The final step is an easy, quick high-gloss polishing with Dura-Polish DIA.

PN 0552



**Dura-Polish DIA** is impregnated with micro-fine diamond particles (67 % by weight) whose structure and distribution are matched to the ceramic fillers of CERAMAGE. This helps to achieve long-lasting gloss and durable composite surfaces in a very short time.

PN 0554



**SHOFU DENTAL GMBH**

Am Brüll 17 · 40878 Ratingen / Germany  
Phone: +49 (0) 21 02 / 86 64 - 0 · Fax: +49 (0) 21 02 / 86 64 - 64  
E-mail: info@shofu.de · www.shofu.de