

Strong bonds without mechanical retention





SHOFU Universal Primer

Easy, efficient and universal

SHOFU Universal Primer reliably bonds resins to all dental alloys, zirconia and alumina


In combination with SHOFU Universal Pre-Opaque, SHOFU Universal Primer creates durable bonds, even in a moist environment, and allows you to apply resins to metal frameworks without any mechanical retention.

For strong bonds

- between light-cured opaque resins and frameworks
- in the fabrication and repair of resin-veneered restorations
- in the intraoral and extraoral repair of PFM restorations
- between denture resins and metal bases
- between resins and orthodontic appliances



SHOFU Universal Primer

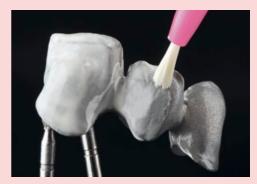

SHOFU Universal Pre-Opaque

Benefits

- Quick and easy handling without any additional equipment
- Colourless, no grey areas at the resin-substrate interface
- Reliable bonding, even without the use of retention beads
- For all dental alloys and zirconia-based frameworks
- For intraoral and extraoral use



SHOFU Universal Primer, the successor to the proven SHOFU MZ Primer Plus, is the universal solution for bonding resins to CAD/CAM metal and zirconia frameworks, thanks to its wide range of indications and quick and easy handling without any further accessories or heat treatments.



SHOFU Universal Opaque

SHOFU Universal Pre-Opaque and SHOFU Universal Primer form a well-matched bonding system, which is extremely resistant to hydrolysis and thermal stress.

After bonding, simply apply SHOFU Universal Opaque, available in 31 shades, in the tooth shade desired.

Order information



SHOFU Universal Primer

Contents 5 ml. PN Y0060



Related products



SHOFU Universal Pre-Opaque

Universal Pre-Opaque is highly flowable and bonds strongly to alloys and Universal Opaque.

Contents: 2 ml, PN 2111



SHOFU Universal Opaque

Light-cured, ready-to-use opaque paste system, including 16 common tooth shades plus Whitening, NCC and Effect shades.

