


For zirconia and lithium disilicate

The three-step diamond impregnated polishing system smoothens the restorations, removes surface scratches made during contouring and provides a high-gloss polish. ZiLMaster is available in a variety of shapes and shanks for both extra and intraoral application.

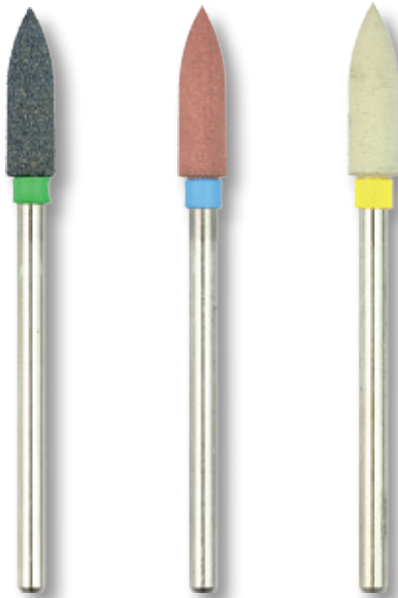

## Indications

Systematic polishing of zirconia and lithium disilicate restorations:

- Finishing/Prepolishing: ZiLMaster Coarse with green banded shank
- Polishing: ZiLMaster Medium with blue banded shank
- Super-polishing: ZiLMaster Fine with yellow banded shank


## highlights

- High gloss in three polishing steps
- Optimal ratio of speed and contact pressure
- Precise concentricity
- Low heat generation
- No polishing paste is needed

PN	Quantity	Product
INDIVIDUAL ARTICLES		
0640-0643, 0646-0649, 0652-0655, 0658-0663, 0665-0668, 0670-0673	3 pcs.	ZiLMaster Coarse, Medium, Fine shapes: bullet, knife, cup, minipoint, KN7, WH6, HP and CA shank
0645, 0651, 0657	1 pc.	ZiLMaster Coarse, Medium, Fine shape: unmounted wheel + HP mandrel
ASSORTMENTS		
0644	1 box	ZiLMaster Coarse, HP shank, 1 each bullet, minipoint, KN7, WH6
0650	1 box	ZiLMaster Medium, HP shank, 1 each bullet, minipoint, KN7, WH6
0656	1 box	ZiLMaster Fine, HP shank, 1 each bullet, minipoint, KN7, WH6
0664	1 box	ZiLMaster Coarse, CA shank, 1 each bullet, knife, cup, minipoint
0669	1 box	ZiLMaster Medium, CA shank, 1 each bullet, knife, cup, minipoint
0674	1 box	ZiLMaster Fine, CA shank, 1 each bullet, knife, cup, minipoint

**Polishers:**  
**Min<sup>-1</sup>** Max. 20,000

**Unmounted wheel:**  
**Min<sup>-1</sup>** Max. 15,000

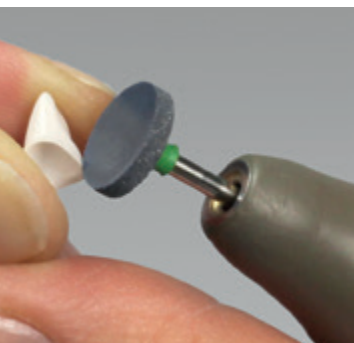




	Shape	Bullet	Knife	Cup	Minipoint	KN7	WH6	Unmounted wheel
	ISO No.	040	090	065	030	125	120	170
ZiLMaster Coarse	HP 3 pcs.	0640			0641	0642	0643	
	HP 1 pc.							0645
	CA 3 pcs.	0658	0662	0663	0661			
ZiLMaster Medium	HP 3 pcs.	0646			0647	0648	0649	
	HP 1 pc.							0651
	CA 3 pcs.	0665	0659	0668	0666			
ZiLMaster Fine	HP 3 pcs.	0652			0653	0654	0655	
	HP 1 pc.							0657
	CA 3 pcs.	0670	0672	0673	0671			

## ZiLMaster - Adjustment kits

The complete system for contouring, finishing and polishing zirconia and lithium disilicate restorations is available in HP and CA shanks. Dura-Green DIA contours the restorations, the three-step system ZiLMaster smoothens the restorations and provides a high gloss.

### HP shank - contents

Dura-Green DIA: 1 each KN7, TC4, WH6  
 ZiLMaster Coarse: 1 each bullet, KN7, WH6  
 ZiLMaster Medium: 1 each bullet, KN7, WH6  
 ZiLMaster Fine: 1 each bullet, KN7, WH6

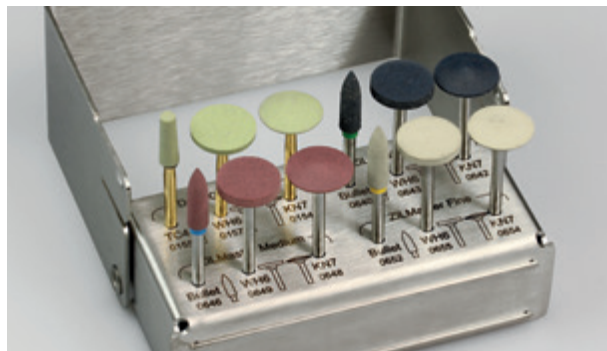

### CA shank - contents

Dura-Green DIA: 1 each FL3, RD3  
 ZiLMaster Coarse: 1 each minipoint, knife, cup  
 ZiLMaster Medium: 1 each minipoint, knife, cup  
 ZiLMaster Fine: 1 each minipoint, knife, cup



PN	Quantity	Product
0680	1 pc.	ZiLMaster Adjustment Kit, HP shank
0681	1 pc.	ZiLMaster Adjustment Kit, CA shank